

CDRPanX Pre-Installation Worksheet


CDR PanX Panoramic System from Schick Technologies, Inc.

Instructions for PanX Installers


■ CDRPANX SYSTEM AND WORKSTATION

CDRPanX: 115v, 20a outlet located within 5 ft (1.5 m) of either side of column centerline.
 Computer: 115v, 15a outlet located within 100 ft (30 m) of CDRPanX.

Top View

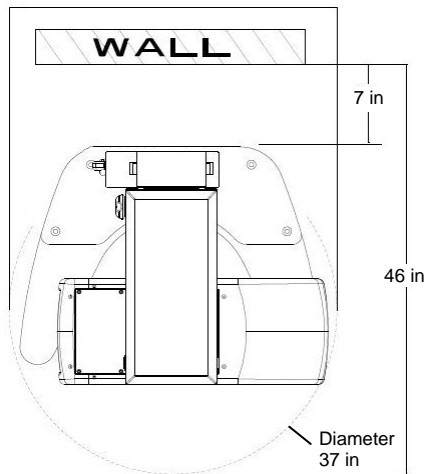


Front View

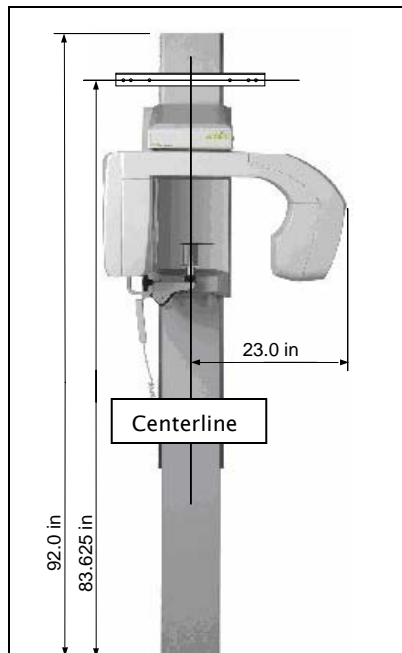

Wall Bracket

Details of the wall bracket are shown below. All items shown on this page are supplied for reference, but are not drawn to scale.

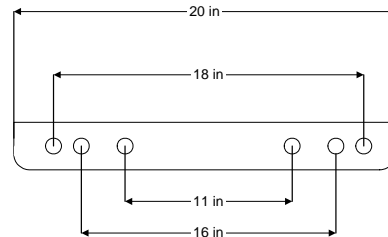

Note for Metal-Stud Sheetrock Walls

For these installations, cabinet hanging screws, similar to #8 fine-threaded self-tapping zinc sheet metal screws, are appropriate. Use a minimum of six screws to install a wood header and then fasten lag bolts from the PanX unit into the wall holder.

IMPORTANT! The above recommendation is intended only as guidance. It is the installer's responsibility to ensure that the wall where the PanX will be installed is sufficient to support and operate the system safely. Refer to the PanX Installation Guide for additional information.